

Europäisches Patentamt
European Patent Office
Office européen des brevets



(11) EP 0 928 832 A2

(12) EUROPEAN PATENT APPLICATION

(43) Date of publication:
14.07.1999 Bulletin 1999/28

(51) Int Cl. 6: C12Q 1/68

(21) Application number: 99300170.0

(22) Date of filing: 11.01.1999

(84) Designated Contracting States:
AT BE CH CY DE DK ES FI FR GB GR IE IT LI LU
MC NL PT SE
Designated Extension States:
AL LT LV MK RO SI

• Kelly, Stephen James
Northwich, Cheshire CW9 7RA (GB)
• Robertson, Nancy Hastings
Northwich, Cheshire CW9 7RA (GB)

(30) Priority: 13.01.1998 GB 9800536

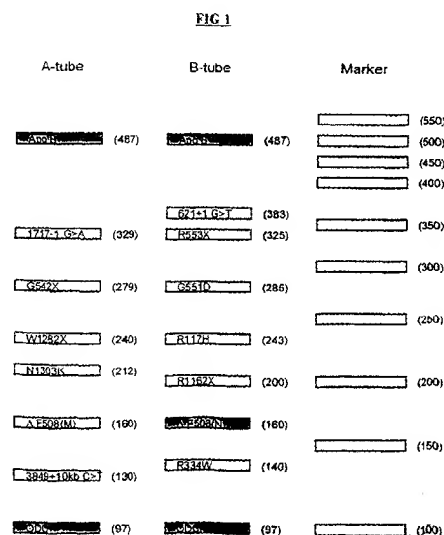
(71) Applicant: ZENECA LIMITED
London W1Y 6LN (GB)

(74) Representative:
Phillips, Neil Godfrey Alasdair et al
Intellectual Property Department
ZENECA Pharmaceuticals
Mersey
Alderley Park
Macclesfield Cheshire SK10 4TG (GB)

(72) Inventors:
• Weston, Susan Louise
Northwich, Cheshire CW9 7RA (GB)

(54) Cystic fibrosis test based on the detection of mutations in the CFRE gene by ARMS

(57) A method for detecting the presence or absence of twelve mutations in the cystic fibrosis transmembrane conductor regulator (CFTR) gene, which method comprises contacting sample genomic DNA from an individual in two separate reaction vessels with allele specific primer sets for (A) 1717-1 G>A, G542X, W1282X, N1303K, AF508(M), 3849+10kb C>T mutations and (B) the 621+1 G>T, R553X, G551D, R117H, R1162X and R334W mutations respectively, in the presence of appropriate nucleotide triphosphates and an agent for polymerisation, such that each diagnostic primer is extended only when the relevant mutation is present in the sample; and detecting the presence or absence of CFTR gene alleles by reference to the presence or absence of diagnostic primer extension product (s).



EP 0 928 832 A2

Description

[0001] Cystic Fibrosis (CF) is the most common fatal autosomal recessive disease affecting Caucasian populations. CF has an incidence of 1 in 2 000 to 3 000 births depending on population group¹ and this indicates a carrier frequency of around 1 in 25 (i.e. 4% of the population). The prognosis for an affected child with CF is a median life expectancy currently estimated to be 40 years². Since the discovery of the CFTR gene³ over 600 CF associated mutations have been identified. The majority of these are presumed to be disease-causing but most are rare so many laboratories test for a limited number with an emphasis on those predominating in their local populations.

[0002] Detection methods used to date have included DNA sequencing, DNA enzyme immunoassay (Sanguinol et al, Int.J.Clin.Lab.Res., 1995, 25, 142-145), multiplex DGGE analysis (Costes et al, Hum.Mutat., 1993, 2, 185-191), and the use of the polymerase chain reaction (PCR) in conjunction with allele-specific oligonucleotide probes (PCR-ASO).

[0003] In addition, ARMS is an established technique that enhances the utility of the PCR for the detection of mutations and other polymorphisms in DNA⁴ (see our European Patent No. 0 332 435). ARMS also provides the means to distinguish between homozygotes and heterozygotes for any given allelic variation. The principle of the method is that under appropriate conditions an oligonucleotide which is not matched to target (genomic) DNA sequence at its 3'-end will not be extended by *Taq* DNA polymerase. The 3'-end of an ARMS primer therefore confers its allele-specificity. Hence, an ARMS product is only generated if the primer is complementary to its target at the 3'-end under the appropriate conditions. We⁵ and others^{6,7} have described the use of ARMS to test for Δ F508, the most common mutation of the CFTR gene³. We have also shown that four ARMS analyses can be performed simultaneously by multiplexing allele-specific primers⁸ and our observations have been confirmed by Fortina *et al*⁹.

[0004] However the need still exists for further CF tests, in particular for rapid and reliable methods for the extensive investigations routinely performed in many laboratories. We have now devised and validated a two-tube multiplex ARMS test, hereinafter referred to as the CF(12)m test, which detects 12 of the most prevalent CF mutations¹⁰ simultaneously. The CFTR gene mutations that are detected by the test are 1717-1G>A, G542X, W1282X, N1303K, Δ F508, 3849+10kb C>T, 621+1 G>T, R553X, G551D, R117H, R1162X and R334W which are described in reference 10 and papers cited therein. The test also distinguishes between CF Δ F508 heterozygotes and homozygotes. The CF Genetic Analysis Consortium data¹⁰ allows the minimum detection capability of the test to be calculated by country, as shown in Table 1. Similarly calculated continental values are shown in Table 2. A further observation from these data is the very high minimum detection capability for Ashkenazic Jews which is calculated to be 95.4%.

[0005] Therefore in a first aspect of the present invention we provide a method for detecting the presence or absence of twelve mutations in the cystic fibrosis transmembrane conductor regulator (CFTR) gene, which method comprises contacting sample genomic DNA from an individual in two separate reaction vessels with allele specific primer sets for (A) 1717-1 G>A, G542X, W1282X, N1303K, Δ F508(M), 3849+10kb C>T mutations and (B) the 621+1 G>T, R553X, G551D, R117H, R1162X and R334W mutations respectively, in the presence of appropriate nucleotide triphosphates and an agent for polymerisation, such that each diagnostic primer is extended only when the relevant mutation is present in the sample; and detecting the presence or absence of CFTR gene alleles by reference to the presence or absence of diagnostic primer extension product(s).

[0006] Allele specific primer sets A and B as defined above represent further independent aspects of the invention.

[0007] Advantageously, primer set B also includes an allele specific primer for the normal Δ F508 allele (Δ F508(N)), i.e. DNA sequence unaffected by the Δ F508 mutation.

[0008] The test sample of nucleic acid is conveniently extracted from a sample of blood, mouthwash or other body fluid from an individual. Also buccal-swab samples may be used. It will be appreciated that the test sample may equally be a nucleic acid sequence corresponding to the sequence in the test sample. That is to say that all or a part of the region in the sample nucleic acid may firstly be amplified using any convenient technique such as PCR before use in the method of the invention. The methods of the invention are conveniently effected using about 20-100 nanograms of template DNA.

[0009] Any convenient enzyme for polymerisation may be used provided that it does not affect the ability of the DNA polymerase to discriminate between normal and mutant template sequences to any significant extent. Examples of convenient enzymes include thermostable enzymes which have no significant 3'-5' exonuclease activity, for example *Taq* DNA polymerase, particularly "Ampli *Taq* Gold"TM DNA polymerase (PE Applied Biosystems), Stoffel fragment, or other appropriately N-terminal deleted modifications of *Taq* or *Tth* (*Thermus thermophilus*) DNA polymerases.

[0010] We have now devised allele specific primers for the above CFTR gene mutations which have been shown, in primer sets A and B, to detect the specific mutations reliably and robustly. Therefore in further aspects of the invention we provide individual diagnostic primers and sets of diagnostic primers A and B as disclosed below:

Primer Set A	Allele specific	Tube	DNA sequence (5' to 3')
W1282X Forward	Y	A	TCTTGGGATTCAATAACTTTGCAACAGTCA
1717-1G>A Forward	Y	A	TACTAAAAGTGACTCTCTAATTTTCTATTTTGGT AATTA
G542X Forward	Y	A	AGTTTGCAGAGAAAGACAATATAGTTCTCT
N1303K Reverse	Y	A	TGATCACTCCACTGTTTCATAGGGATCCATC
ÆF508 Mutant Reverse	Y	A	GTATCTATATTCATCATAGGAAACACCATT
3849+10kb C>T Reverse	Y	A	ACATTTCTTTTCAGGGTGTCTGACTAA

Primer Set B	Allele specific	Tube	DNA sequence (5' to 3')
ÆF508 Wild type Reverse	Y	B	GTATCTATATTCATCATAGGAAACACCACA
621+1 G>T Reverse	Y	B	TGCCATGGGGCCTGTGCAAGGAAGTATTGA
R117H Reverse	Y	B	AGCCTATGCCCTAGATAAATCGCGATAGACT
R334W Forward	Y	B	CCTATGCACTAATCAAAGGAATCATCCTGT
G551D Reverse	Y	B	GCTAAAGAAATTCTTGCTCGTTGTT
R553X Reverse	Y	B	GACTGACTGACTGACTGACTCTGACTGACTTATT CACCTTGCTAAAGAAATTCTTGCTGA
R1162X Forward	Y	B	TATTTTATTTTCAGATGCGATCTGTGAGTT

[0011] Although the above primers have been optimised for use in the diagnostic methods of the invention, it will be appreciated that some modification of some or all of the above primers may be possible without adversely affecting their performance in the methods of the invention. Such modifications may be determined by the molecular biologist of ordinary skill in comparative experiments. Whilst we do not wish to be bound by theoretical considerations, in general, one or more of the nucleotides at the 5' end of a primer, may be substituted for other (non-complementary) nucleotides, for example about 1 in 4, 1 in 5, 1 in 6, 1 in 7, such as 1 in 8, 1 in 9 or 1 in 10 of the nucleotides may be substituted. Nearer to the 3' end of the primer the scope for modification is more limited, within about six bases of the 3' end, one or possibly two nucleotides may be substituted for other (non-complementary) nucleotides. The 3' terminal nucleotide of the primer should not be modified.

[0012] Additions to and deletions from the above primer sequences are also possible but not preferred. Again a molecular biologist may perform comparative experiments to reveal what modifications are practicable within the scope of the present invention.

[0013] The primers may be manufactured using any convenient method of synthesis. Examples of such methods may be found in standard textbooks, for example "Protocols For Oligonucleotides And Analogues: Synthesis And Properties," Methods In Molecular Biology Series; Volume 20; Ed. Sudhir Agrawal, Humana ISBN: 0-89603-247-7; 1993; 1st Edition.

[0014] In many situations, it will be convenient to use the sets of diagnostics primer of the invention in combination with further amplification primers in one or more cycles of PCR amplification. A convenient example of this aspect is set out in our European patent number EP-B1-0332435. The further (amplification) primer is either a forward or a reverse common primer. Hot-start PCR is conveniently employed.

[0015] Preferred individual primers and sets of further primers are set out below.

EP 0 928 832 A2

Primer	Allele specific	Tube	DNA sequence (5' to 3')
W1282X Reverse	N	A	GAATTCCTCAAACTTTAGAGACATC
1717-1G>A/G542X Reverse	N	A	TAATCTCTACCAAATCTGGATACTATACC
N1303K Forward	N	A	AATTTGAGAGAACTTGATGGTAAGTACA
ÆF508 Mutant Forward	N	A	CCAGACTTCACTTCTAATGATGATTATGGG
3849+10kb C>T Forward	N	A	TTGTGGATCAAATTCAGTTGACTTGTCATC
ÆF508 WT Forward	N	B	GACTTCACCTCTAATGATGATTATGGGAGA
621+1 G>T/R117H Forward	N	B	GTTTCACATAGTGTATGACCCTCTATATACACTCATT
R334W Reverse	N	B	TTTGTATTGCTCCAAGAGAGTCATACCA
G551D/R553X Forward	N	B	TAAAATTGGAGCAATGTTGTTTTGACC
R1162X Reverse	N	B	TTTGTCTGTGAGATCTTTGACAGTCATT

[0016] In a preferred aspect of the invention primer mixes A and B are provided which mixes A and B comprise both allele specific and amplification primers as follows:

Primer Mix A	Allele specific	Tube	DNA sequence (5' to 3')
W1282X Forward	Y	A	TCTTGGGATTCAATAACTTTGCAACAGTCA
W1282X Reverse	N	A	GAATTCCTCAAACTTTAGAGACATC
1717-1G>A Forward	Y	A	TACTAAAAGTGACTCTCTAATTTTCTATTTTGGT AATTA
G542X Forward	Y	A	AGTTTGACAGAGAAAGACAATATAGTTCTCT
1717-1G>A/G542X Reverse	N	A	TAATCTCTACCAAATCTGGATACTATACC
N1303K Reverse	Y	A	TGATCACTCCACTGTTTCATAGGGATCCATC
N1303K Forward	N	A	AATTTGAGAGAACTTGATGGTAAGTACA
ÆF508 Reverse	Y	A	GTATCTATATTCATCATAGGAAACACCATT
ÆF508 Forward	N	A	CCAGACTTCACTTCTAATGATGATTATGGG
3849+10kb C>T Reverse	Y	A	ACATTTCTTTTCAAGGGTGTCTGACTAA
3849+10kb C>T Forward	N	A	TTGTGGATCAAATTCAGTTGACTTGTCATC

Primer Mix B	Allele specific	Tube	DNA sequence (5' to 3')
ÆF508 Wild type Reverse	Y	B	GTATCTATATTCATCATAGGAAACACCACA
ÆF508 Forward	N	B	GACTTCACCTCTAATGATGATTATGGGAGA
621+1 G>T Reverse	Y	B	TGCCATGGGGCCTGTGCAAGGAAGTATTGA
R117H Reverse	Y	B	AGCCTATGCCTAGATAAATCGCGATAGACT
621+1 G>T/R117H Forward	N	B	GTTTCACATAGTGTATGACCCTCTATATACACTC ATT
R334W Forward	Y	B	CCTATGCACTAATCAAAGGAATCATCCTGT
R334W Reverse	N	B	TTTGTATTGCTCCAAGAGAGTCATACCA
G551D Reverse	Y	B	GCTAAAGAAATTCCTGCTCGTTGTT
R553X Reverse	Y	B	GACTGACTGACTGACTGACTCTGACTGACTTATT CACCTTGCTAAAGAAATTCCTGCTGA
G551D/R553X Forward	N	B	TAAAATTGGAGCAATGTTGTTTTGACC
R1162X Forward	Y	B	TATTTTATTTTCAAGATCGATCTGTGAGTT
R1162X Reverse	N	B	TTTGTCTGTGAGATCTTTGACAGTCATT

[0017] Any convenient control primer(s) may be used to confirm correct functioning of the PCR. Two independent control reactions may be provided with each primer set. For example, one amplifies a region of the apolipoprotein (Apo B) gene, the other amplifies a region of the ornithine decarboxylase (ODC) gene. The Apo B reaction is designed to amplify less efficiently than the ARMS reactions ensuring that this control reaction will fail preferentially if PCR variables are not within the defined tolerances of the test. The Apo B control therefore safeguards against a false negative result

if the test is performed under sub-optimal conditions. Both control product bands are preferably present for a reliable diagnosis to be made.

[0018] Preferred control primers include one or more of the following:

Apo B control Forward	N	A & B	GAGCACAGTACGAAAAACCACT
Apo B control Reverse	N	A & B	AAACTTTTACAGGGATGGAGAACG
ODC control Forward	N	A & B	AGAGGATTATCTATGCAAATCCTTGTAACC
ODC control Reverse	N	A & B	TCAACTTCACTATCAAAAGTCATCATCTAG

[0019] It will be appreciated that any of the primers set out hereinbefore may be modified as outlined in respect of primer sets A and B.

[0020] A variety of methods may be used to detect the presence or absence of diagnostic primer extension products and/or amplification products. These will be apparent to the person skilled in the art of nucleic acid detection procedures. Preferred methods avoid the need for radiolabelled reagents. Particular detection methods include those which provide size or sequence differentiation, preferably gel electrophoresis based methods. For example, primer extension products are separated by agarose gel electrophoresis and visualised by ethidium bromide staining/uv transillumination. Individual mutations are conveniently identified by size comparison of the products against an appropriate DNA ladder marker. Alternative detection methods include capillary zone electrophoresis (CZE).

[0021] One or more of the diagnostic primer sets/mixes of the invention may be conveniently packaged with instructions for use in the method of the invention and appropriate packaging and sold as a kit. The kits will conveniently include one or more of the following: corresponding amplification primers, appropriate nucleotide triphosphates, for example dATP, dCTP, dGTP, dTTP, a suitable polymerase as previously described, and optimised buffer solution.

[0022] A preferred kit of the invention comprises primer mixes A and B (as defined above) in separate containers, taq DNA polymerase (preferably Ampliqa Gold), and optimised enzyme dilution buffer.

[0023] The invention will now be illustrated but not limited by reference to the following detailed description, References, Examples, Tables and Figures wherein

Table 1 shows the minimum detection capabilities of the CF(12)m ARMS test; average figures for countries testing >100 chromosomes¹⁰.

Table 2 shows the minimum detection capabilities of the CF(12)m ARMS test: analysis by continent derived from data in reference 10. Empty cells indicate that not all centres screened for each mutation detected by the test; therefore the actual detection capabilities for some mutations may be higher.

Table 3 shows the analysis of the 754 chromosomes tested. "Confirmatory typing as detailed in references cited within reference 10: ASO, Allele-specific oligonucleotide hybridisation.

Table 4 shows primer sequences used in the CF(12)m ARMS test.

[0024] Figure 1 shows a diagrammatic representation of electrophoretic mobilities of ARMS, control and size marker gel bands. The respective ARMS and control bands are labelled, all band sizes (bp) are shown in parentheses. The 250 bp marker band (hatched) is more intense to provide a reference position on the gels. In samples which do not contain any of the mutations detected by the test, only the Apo B and ODC control bands (both tubes) and Δ F508 normal band (B-tube) will be visible. these are shown shaded.

References:

[0025]

- Welsh MJ, Tsui L-C, Boat TF, In Scriver CR, Beaudet AL, Sly WL, Valle (eds), The metabolic and molecular basis of inherited disease. McGraw-Hill, New York: 3799-3876, 1990.
- Elborn JS, Shale DJ, Britton JR. Cystic fibrosis current survival and population estimates to the year 2000. *Thorax* 1991;46:881-885.
- Riordan JR, Rommens JM, Kerem B, Alon N, Rozmahel R, Grzelczak Z, *et al.* Identification of the cystic fibrosis gene: cloning and characterization of complementary DNA. *Science* 1989;245:1066-1073.
- Newton CR, Graham A, Heptinstall LE, Powell SJ, Summers C, Kalsheker N, *et al.* Analysis of any point mutation in DNA. The amplification refractory mutation system (ARMS). *Nucl Acids Res* 1989;17:2503-2516.
- Newton CR, Heptinstall LE, Summers C, Super M, Schwarz M, Graham A, *et al.* Detection of delta-F508 deletion by amplification refractory mutation system. *Lancet* 1990;i:1217-1219.

6. Ballabio A, Gibbs RA, Caskey CT. PCR test for cystic fibrosis deletion. *Nature* 1990;**343**:220.
7. Wagner M, Schloesser M, Reiss J. Direct gene diagnosis of cystic fibrosis by allele-specific polymerase chain reactions. *Mol Biol Med* 1990;**7**:359-364.
8. Ferrie RM, Schwarz MJ, Robertson NH, Vaudin S, Super M, Malone G, *et al*. Development, multiplexing and application of ARMS tests for common mutations in the CFTR gene. *Am J Hum Genet* 1992;**51**:251-62.
9. Fortina P, Conant R, Monokian G, Dotti G, Parrella T, Hitchcock W, *et al*. Non-radioactive detection of the most common mutations in the cystic fibrosis transmembrane conductance regulator gene by multiplex allele-specific polymerase chain reaction. *Hum Genet* 1992;**90**:375-378.
10. Kazazian HH. Population variation of common cystic fibrosis mutations. *Hum Mutat* 1994;**4**:167-177.
11. Kwok S, Higuchi R. Avoiding false positives with PCR. *Nature* 1989;**339**:237-238.
12. Super M, Schwarz MJ, Malone G, Roberts T, Haworth A, Dermody G. Active cascade testing for carriers of cystic fibrosis gene. *BMJ* 1994;**308**:1462-1467.
13. Wald NJ. Couple screening for cystic fibrosis. *Lancet* 1991;**338**:1318-1319.
14. Gilfillan A, Axton R, Brock DJH. Mass screening for cystic fibrosis heterozygotes: Two assay systems compared. *Clin Chem* 1994;**40**:197-199.
15. Livingstone J, Axton FA, Gilfillan A, Mennie M, Compton M, Liston WA, *et al*. Antenatal screening for cystic fibrosis: a trial of the couple model. *BMJ*, 1994; **308**:1459-1462.
16. Findlay I, Cuckle H, Lilford RJ, Rutherford RJ, Quirke P, Lui S. Screening sperm donors for cystic fibrosis. *BMJ* 1995;**310**:1533.
17. Cuppens H, Cassiman JJ. A quality control study of CFTR mutation screening in 40 different European laboratories: The European concerted action on cystic fibrosis. *Eur J Hum Genet* 1995;**3**:235-245.

Example 1

- [0026] For clinical applications it is essential that the results obtained using PCR based tests are both reliable and reproducible. It is generally understood that the efficiency of PCR can be affected by a number of variables and it is important to understand which ones might affect test performance. We have therefore investigated the effects of changes of primer, DNA template and Taq DNA polymerase concentrations and of PCR annealing temperatures. These studies have defined the conditions under which neither false positive nor false negative results are produced even if the test is performed under sub-optimal conditions. We have unequivocally typed DNA samples that were obtained from different laboratories and that were prepared by a variety of methods.

Methods

DNA samples

- [0027] The panel of samples was selected to include normal DNAs, DNA from CF Δ F508 homozygotes and several examples of each mutation for which the test was designed to detect. All DNA samples were prepared from EDTA/blood. DNA samples from external source were prepared and typed independently using standard recognised procedures and typed by the methods outlined in table 3. DNA samples from CF unaffected individuals were prepared as described previously⁸.

Primer design

- [0028] In designing the ARMS primers it was important to ensure that a false result could not arise due to other DNA sequence variations at the same site or in the vicinity of the mutations tested for. For example, Δ F508 and non- Δ F508 alleles should not be confused for either the mutant Δ I507, or the benign F508C alleles. This discrimination was achieved after particular consideration and the appropriate choice of orientation of each primer with respect to the direction of transcription of the CFTR gene. Careful attention was also given to the inclusion of additional base-pair mismatches between each primer and the genomic DNA sequence and also the length and concentration of each primer. In combining primers for multiplex analysis, many primers are included in one reaction mix. It was therefore necessary to minimise any primer interactions that might affect the test performance. The ARMS primer sequences are shown in Table 4.

Test design

- [0029] The CF(12)m test consists of two tubes, the A and the B tube. The A tube contains ARMS primers specific for the 1717-1G>A, G542X, W1282X, N1303K, AF508, and 3849+10kb C>T mutations. The B tube contains ARMS

EP 0 928 832 A2

primers specific for the 621+1 G>T, R553X, G551D, R117H, R1162X and R334W mutations. The B tube also contains an ARMS primer specific for the normal Δ F508 allele. There are also two control reactions in each tube. One amplifies a region of the apolipoprotein (Apo B) gene, the other amplifies a region of the ornithine decarboxylase (ODC) gene. The Apo B reaction is designed to amplify less efficiently than the ARMS reactions ensuring that this control reaction will fail preferentially if PCR variables are not within the defined tolerances of the test. The Apo B control therefore safeguards against a false negative result if the test is performed under sub-optimal conditions. Both control product bands should be present for a reliable diagnosis to be made. The concentrations of the component ARMS primers in each tube were adjusted to amplify their respective targets with equal efficiency. The primer DNA sequences are shown in Table 4. The remaining constituents of each tube comprise ARMS buffer; 10mM Tris-HCl, (pH 8.3) 1.2mM MgCl₂, 50mM KCl, 0.01% gelatin, dNTPs (100 μ M each). Genomic DNA and *Taq* or *Amplitaq Gold* DNA polymerase are added at the start of the test.

[0030] The CF(12)m test output is shown (Figure 1)

Test method

[0031] The test is performed by adding genomic DNA prepared from EDTA/blood to each tube which is then heated at 94°C, 5 minutes. *Taq* DNA polymerase (2.5 units) is added to each tube, whilst maintaining the temperature at 94°C, and thermal cycling (35 cycles of 94°C, 1 minute; 58°C, 2 minutes; 72°C, 1 minute with a final extension at 72°C, 10 minutes) is initiated. Alternatively, both *Amplitaq Gold* DNA polymerase (2.5 units) and genomic DNA prepared from EDTA/blood are added to each tube and thermal cycling (as described above) initiated following a 94°C, 20 minutes incubation during which the DNA polymerase is activated.

[0032] The amplification products are separated by agarose gel electrophoresis against a 50 base-pair DNA ladder marker (Pharmacia Biotech) and visualised by UV transillumination. All procedures generally accepted for avoiding PCR carry-over contamination¹¹ were employed during the development and validation of the test and are recommended when using the test. The presence of the control products and a specific ARMS product, defined by electrophoretic mobility, is diagnostic of the presence of the respective mutant allele and/or normal Δ F508 allele (Figure 1).

Results

[0033] Test performance:

[0034] A series of experiments was performed to determine the tolerance of the CF(12)m ARMS test system to the variation of primer concentrations, annealing temperatures, *Taq* or *Amplitaq Gold* DNA polymerase concentration, the input DNA amount and the method of its extraction. This served to demonstrate the permissible limits of each test variable and to confirm that a 'non-result' or test fail, rather than an incorrect result was observed should any variable fall outside defined tolerances. The CF(12)m test was functional at $\pm 2^\circ\text{C}$ of the standard annealing temperature (58°C) Under sub-optimal conditions (3°C above the standard temperature) all diagnostic product bands were visible but the upper control product band was absent. This feature of the test prevents the occurrence of false negative results due to operation at higher temperatures, since both control product bands must be present to make a diagnosis. The results demonstrate that the test remains functional over a broad range of annealing temperatures. The effects of the primer, *Taq* or *Amplitaq Gold* DNA polymerase and DNA concentrations were assessed in an analogous manner. This defined the respective windows within which both of the control bands and each mutation-specific band was appropriately generated.

[0035] 383 DNA samples prepared from EDTA/blood were analysed using the CF(12)m test, and 377 gave results in agreement with those obtained independently. Six samples (1.57%) were not scored due to the absence of one or both control bands. These were all archival samples known to be of low concentration and degradation of these samples could not be eliminated. Where rare non- Δ F508 compound heterozygotes have been obtained (3849+10kb C>T/W1282X; 3849+10kb C>T/G542X; G542X/N1303K; G542X/W1282X; G551D/R553X; N1303K/1717-1G>A; G542X/1717-1G>A; N1303K/W1282X; R553X/R334W) and analysed, both mutations were correctly identified. Similarly, all Δ F508 compound heterozygotes were correctly typed. We also observed that there was no mistyping when Δ 1507, 1717-2A>G, R1283M, R117C, 3617G/T, 621+2T>C or F508C alleles were present. These data demonstrate that, for the clinical material available to us, one allele does not fail to amplify when combined with another in the same ARMS tube and that one allele does not give rise to artefactual products derived from alternative allele-specific primers. Table 3 provides an analysis of the 754 chromosomes typed using CF(12)m ARMS test.

[0036] Early diagnosis followed by expert management has resulted in an improved prognosis for CF patients. This has led to the evaluation of neonatal screening protocols for CF. There is also a common objective of identifying carriers and ultimately those couples who are at risk of having a child affected with CF that is shared by 'Active cascade'¹², 'Couple'^{13,14} and 'Two step' or 'Stepwise'^{14,15} screening protocols. It has also been highlighted that there are risks associated with *in vitro* fertilisation. Up to ten offspring may be fathered by one sperm donor. If a donor is a CF carrier

there will be a higher than normal chance of at least one of the offspring being affected by the disease. For this reason it was recommended that sperm donors should be routinely tested for mutations in the CFTR gene¹⁶. The primary application of the CF(12)m test is screening individuals who may be carriers of one of the CF mutations 1717-1G>A, G542X, W1282X, N1303K, AF508, 3849+10kb C>T, 621+1 G>T, R553X, G551D, R117H, R1162X and R334W, the most common CF mutations in Caucasians and Ashkenazi Jews. In Europe, CF individuals who are non-AF508 CF compound heterozygotes are rare. As a consequence it has not been possible to analyse for all mutations in combination in the same test. However, nine compound heterozygotes have been sourced and analysed. In each of these cases and in all AF508 CF compound heterozygotes both mutations were correctly identified clearly demonstrating the diagnostic value of the test.

[0037] There is a clear difficulty in distinguishing particular alleles in the routine genetic screening laboratory¹⁷. Cuppens and Cassiman demonstrated that 12.5% of laboratories mistyped the F508C polymorphism as a true mutation and that 12.5% confused the Δ I507 mutation for Δ F508¹⁷. One conclusion of that study was that the accuracy of CFTR typing should be improved¹⁷. The CF(12)m test was designed such that a false result does not arise because of other DNA sequence variations in the vicinity of the mutations tested for. This was achieved through the design and choice of transcriptional orientation of each respective ARMS primer. The CF(12)m test accurately discriminates between Δ F508 and non- Δ F508 alleles which are in turn distinguished from the mutant Δ I507 and the benign F508C alleles. This addresses the problems associated with typing particular alleles identified by Cuppens and Cassiman. Furthermore, findings from another study¹⁴, that compared multiplex ARMS screening for the Δ F508, G551D, G542X and 621+1 G>T alleles⁸ with alternative routine procedures for the same alleles were that multiplex ARMS was the preferred method.

[0038] Standardisation of all test procedures and implementation of appropriate quality control measures can control for some test variables and it is important to understand which ones are likely to affect test performance. The annealing temperature is one key parameter which can affect the performance of any PCR or ARMS reaction. If the reaction temperature falls below the required annealing temperature it is likely that reaction specificity will be affected. Similarly, if the annealing temperature is higher than that required the efficiency of the ARMS reaction may be reduced resulting in PCR failure. There is expected to be some degree of variability between the temperature accuracy and uniformity exhibited by different thermal cyclers. We therefore established the temperature range within which all the component reactions of the test would function (at least $\pm 2^\circ\text{C}$). All thermal cyclers currently available commercially are claimed by the requisite manufacturer to have a temperature accuracy of $\pm 1^\circ\text{C}$ or better. We have similarly determined the range within which other parameters can vary without adversely affecting test function. In defining these we have established the tolerances within which the kit must be assembled. The two key variables that can be influenced by the kit user are the annealing temperature of the ARMS reaction and the amount of the *Taq* DNA polymerase added to individual reaction mixes. We have addressed the function of the test when the upper and lower limits of these variables are applied in concert. From these studies, even with low enzyme and high annealing temperature and vice versa, the control amplicons are produced and the test does not give false results using our panel of DNAs. We therefore conclude that the CF(12)m test is reliable and robust and is unaffected by external variations so long as the kit instructions and recommendations are adhered to. If a reaction should fail for any other reason that is not addressed by the parameters investigated herein, a misdiagnosis would still not occur since a control amplicon would be eliminated prior to any ARMS amplicon. This feature emphasises the importance of the internal controls and their role in avoiding misdiagnoses. There were six archival DNA samples that failed to produce one or both control bands, these may have been degraded and they were not scored. Because the CF(12)m test is primarily a screening test, by implication DNA used in the test will be freshly prepared. The age and degradation of a sample is therefore not believed to be a problem during the test's routine application.

Example 2

Clinical evaluation of the CF(12)m cystic fibrosis DNA diagnostic kit

[0039] Forty DNA samples from thirty nine previously genotyped cystic fibrosis (CF) patients or carriers and one unaffected individual were evaluated independently in a blind study by two genetic testing laboratories involved in CF molecular diagnosis performing 500 diagnoses per year combined. The CF(12)m kit uses multiplexed amplification refractory mutation system (ARMSTM) technology which allows the simultaneous identification of the more prevalent CFTR gene mutations⁵ 1717-1G>A, G542X, W1282X, N1303K, Δ F508, 3849+10kb C>T, 621+1 G>T, R553X, G551D, R117H, R1162X and R334W in one working day. The kit was found to be accurate and reliable.

Methods

Samples

5 [0040] Forty DNA samples from CF patients or carriers and an unaffected individual were typed. The reference method for DNA typing was the analysis of the 27 exons and the intron-exon boundaries of the CFTR gene as described previously⁶. Six further samples bearing the Δ I507, R117C, R347H, D1152H and R117P mutations and one bearing the 1540A->G polymorphism, not detectable with the CF(12)m kit, were also tested to evaluate cross reactivity. Whole blood had been stored at -20°C prior to CF(12)m kit analysis.

DNA analysis using CF(12)m kits

10 [0041] DNA was extracted from thawed blood and analysed. After DNA amplification the reaction mixtures were separated by agarose gel electrophoresis to reveal the diagnostic and amplification control DNA fragments³.

Results

[0042] The extraction of blood samples could be completed in half a day. DNA was successfully amplified from all samples, however, one remained dark brown after washing and could not be amplified initially. This DNA was diluted ten fold after which it amplified successfully.

20 [0043] All samples were analysed in both laboratories. In one lab one sample required a repeat analysis because of the absence of a control DNA fragment. One sample, different in each lab (see Table 1), was not interpreted even after reamplification. In lab 1 sample six had an allele carrying the 1717-1G->A mutation which was amplified only at a low level. In lab 2 sample 26 showed no amplification control band and no mutant ARMS products.

25 [0044] As expected, the mutations 2789+5 G->A, 2176insC, 1078delT and S1235R and the polymorphisms 875+40A/G and 2694A/C identified by the reference method were not detected by the kit, neither were they erroneously typed. This was also found for the Δ I507, R117C, R347H, D1152H, R117P 1540A->G samples [data not shown]. None of the CF mutations detectable by the kit were incorrectly identified as other mutations.

[0045] In all cases with the exception of samples six and 26 where one or more of the 12 mutations detected by the kit were present, the kit accurately identified them. When compared to the reference method the kit detected the same mutations in 34 of the 40 samples (85%), (see Table 1). Four samples carried a rarer mutation that is not detected by the kit. Two other samples were from patients homozygous for the W1282X and G542X mutations. Figure 1 shows a result typical of those produced using the CF(12)m kit.

35 [0046] The CF(12)m kit, uses multiplex ARMS technology which is extremely rapid and allows the screening of at least ten samples for 12 CFTR mutations in a working day. A major feature of the kit therefore is the speed with which the results are delivered. This allows the laboratory staff to undertake additional activities and it increases sample throughput. We found the kit to be reliable and accurate after evaluating it in a blind study in two independent laboratories both of which are experienced in CF molecular diagnostics. The results obtained were in accordance with those described by the manufacturer and demonstrated a first time amplification success rate of 97%.

40 [0047] In one laboratory, the absence or low level of the amplification control DNA fragment was observed and this was improved by the analysis of small batches of samples. This probably reflects reduced denaturation of the DNA polymerase during the hot-start prior to thermal cycling. To achieve consistently reliable results we found that strict adherence to the protocols was required. The control DNA fragments should be visible routinely to avoid the risk of getting a false negative result. The absence of the upper control fragment in conjunction with no mutant allele amplicons indicate that a repeat test is required.

45 [0048] The 12 mutations that the CF(12)m kit can detect are among the most frequent observed in European populations. The diagnostic power of the test, depending on the population under study has been calculated. In the French population generally, 78.8% of CF chromosomes should be detected using the CF(12)m kit, in our hands 85% were detected in both laboratories. This seems satisfactory as screening for very rare mutations does not significantly increase the detection frequency of CF chromosomes.

50 [0049] A potential limitation of the CF(12)m kit is that it does not distinguish between homozygotes and heterozygotes for mutations other than AF508. Although homozygotes for the rarer mutations are usually encountered infrequently, this characteristic of the kit must be considered when interpreting the test result. This was demonstrated in this investigation when the genotype of the W1282X and G542X homozygous patients could not be determined since the kit does not differentiate between one or two mutated alleles except for the AF508 mutation.

55 [0050] The identification of mutations in the CFTR gene is required primarily in patients suspected of suffering from CF and is important for genetic counselling in the patient family. The CF(12)m kit is useful in this case for first line screening and can be complemented by other methods if two mutations are not detected initially. The detection rate

EP 0 928 832 A2

of the kit allows the calculation by Bayesian analysis of the posterior risk of an 'at risk individual', with a negative result. Similar assessments of carrier risk estimation can be calculated for individuals with a family history of CF where affected or obligate heterozygous individuals are not available for mutation typing.

[0051] The CF(12)m kit will also be useful in CF carrier screening programs and of particular use in CF neonatal screening in association with immunoreactive trypsin (IRT) (Wilcken et al, J.Pediat., 1995, 127, 965-970) or pancreatitis-associated protein (PAP) (CR.Acad.Sci.Paris/Life Sciences, 1994, 317, 561-564).

Table 1

Country	CF chromosomes (%)
Albania	75.0
Austria	72.2
Belgium	84.5
Bulgaria	67.1
Czech Republic	80.6
Denmark	91.0
France	78.8
Germany	78.6
Hungary	73.0
Greece	68.7
Ireland	79.1
Italy	60.3
Netherlands	83.5
Poland	65.0
Portugal	48.4
Russia	53.4
Slovenia	69.1
Spain	63.2
Sweden	56.0
Switzerland	80.3
United Kingdom	81.3

Table 2

	Total CF chromo- somes screened	R117H	621+1G> T	R334 W	AEF508	1717-1G>A	G542X
Europe, North	21,154	62	97	18	14,866	160	439
Europe, South	7,281	3	37	21	4,007	65	259
America, North	10,438	61	154	12	6,900	44	234
America, South/ Central	758				342		38
Australasia	3,095	7	27	2	2,309	12	56
Asia, mainly Middle East	608	0			173	3	27
Africa	515				351		9
Totals	43,849	133	315	53	28,948	284	1,062

	G551D	R553X	R1162X	3849+10 kbC>T	W1282X	N1303K	Minimum Detection Capability (%)
Europe, North	356	165	36	23	120	209	78.2
Europe, South	37	44	68	8	43	179	65.5
America, North	206	96	19	57	245	130	78.2
America, South/ Central	1	5				11	52.4
Australasia	117	11	?		6	23	83.1
Asia, mainly Middle East	0	0		16	120	29	60.5
Africa		1			2	8	72.0
Totals	717	322	125	104	536	589	75.7

Table 3

Mutation	Independent typing method ^a	Totals
1717-1G>A	ASO	16
G542X	ASO	10
W1282X	ASO	16
N1303K	ASO	12
AEF508	Electrophoresis	89
3849+10kb C>T	Digest (HphI)	11
621+1 G>T	Digest (MseI)	7

EP 0 928 832 A2

Table 3 (continued)

Mutation	Independent typing method ^a	Totals
R553X	Digest (HincII)	15
G551D	Digest (NdeI)	16
R117H	ASO	13
R1162X	Digest (DdeI)	11
R334W	Digest (MspI)	6
Other/none		532
Number of samples		377
Total number of chromosomes		754

Table 4

Primer	Allele specific	Tube	DNA sequence (5' to 3')
Apo B control Forward	N	A & B	GAGCACAGTACGAAAAACCACCT
Apo B control Reverse	N	A & B	AAACTTTTACAGGGATGGAGAACG
ODC control Forward	N	A & B	AGAGGATTATCTATGCAAATCCTTGT AACC
ODC control Reverse	N	A & B	TCAACTTCACTATCAAAAGTCATCAT CTAG
W1282X Forward	Y	A	TCTTGGGATTCAATAACTTTGCAACA GTCA
W1282X Reverse	N	A	GAATTCCTAACTTTTAGAGACATC
1717-1G>A Forward	Y	A	TACTAAAAGTGACTCTCTAATTTTCTA TTTTTGGTAATTA
G542X Forward	Y	A	AGTTTGCAGAGAAAGACAATATAGTT CTCT
1717-1G>A/G542X Reverse	N	A	TAATCTCTACCAAATCTGGATACTAT ACC
N1303K Reverse	Y	A	TGATCACTCCACTGTTTCATAGGGATC CATC
N1303K Forward	N	A	AATTTGAGAGAACTTGATGGTAAGTA CA
ÆF508 Reverse	Y	A	GTATCTATATTCATCATAGGAAACAC CATT
ÆF508 Forward	N	A	CCAGACTTCACTTCTAATGATGATTA TGGG
3849+10kb C>T Reverse	Y	A	ACATTTCTTTTCAGGGTGTCTGACTA A
3849+10kb C>T Forward	N	A	TTGTGGATCAAATTTCAAGTTGACTTG TCATC
ÆF508 Wild type Reverse	Y	B	GTATCTATATTCATCATAGGAAACAC CACA
ÆF508 Forward	N	B	GACTTCACTTCTAATGATGATTATGG GAGA
621+1 G>T Reverse	Y	B	TGCCATGGGGCCTGTGCAAGGAAGT ATTGA
R117H Reverse	Y	B	AGCCTATGCCTAGATAAATCGCGATA GACT
621+1 G>T/R117H Forward	N	B	GTTTCACATAGTGATGACCCTCTAT ATACACTCATT
R334W Forward	Y	B	CCTATGCACTAATCAAAGGAATCATC CTGT
R334W Reverse	N	B	TTTGTTTATTGCTCCAAGAGAGTCAT ACCA
G551D Reverse	Y	B	GCTAAAGAAATTCTTGCTCGTTGTT
R553X Reverse	Y	B	GACTGACTGACTGACTGACTCTGAC TGACTTATTCACCTTGCTAAAGAAAT TCTTGCTGA
G551D/R553X Forward	N	B	TAAAATTGGAGCAATGTTGTTTTGA CC
R1162X Forward	Y	B	TATTTTATTTTCAAGATGCGATCTGTG AGTT
R1162X Reverse	N	B	TTTGTCTGTGAGATCTTTGACAGTCA TTT

Annex to the description

[0052]

5

SEQUENCE LISTING

10

(1) GENERAL INFORMATION:

(i) APPLICANT:

15

- (A) NAME: Zeneca Ltd
- (B) STREET: 15 Stanhope Gate
- (C) CITY: London
- (D) STATE: Greater London
- (E) COUNTRY: England
- (F) POSTAL CODE (ZIP): W1Y 6LN
- (G) TELEPHONE: 0171 304 5000
- (H) TELEFAX: 0171 304 5151
- (I) TELEX: 0171 834 2042

20

(ii) TITLE OF INVENTION: SEQUENCES

25

(iii) NUMBER OF SEQUENCES: 27

(iv) COMPUTER READABLE FORM:

30

- (A) MEDIUM TYPE: Floppy disk
- (B) COMPUTER: IBM PC compatible
- (C) OPERATING SYSTEM: PC-DOS/MS-DOS
- (D) SOFTWARE: PatentIn Release #1.0, Version #1.30 (EPO)

35

(vi) PRIOR APPLICATION DATA:

- (A) APPLICATION NUMBER: GB 9800536.6
- (B) FILING DATE: 13-JAN-1998

40

(2) INFORMATION FOR SEQ ID NO: 1:

(i) SEQUENCE CHARACTERISTICS:

45

- (A) LENGTH: 23 base pairs
- (B) TYPE: nucleic acid
- (C) STRANDEDNESS: single
- (D) TOPOLOGY: linear

50

(ii) MOLECULE TYPE: other nucleic acid

55

(xi) SEQUENCE DESCRIPTION: SEQ ID NO: 1:

GAGCACAGTA CGAAAACCA CCT

23

(2) INFORMATION FOR SEQ ID NO: 2:

- 5 (i) SEQUENCE CHARACTERISTICS:
 (A) LENGTH: 24 base pairs
 (B) TYPE: nucleic acid
 (C) STRANDEDNESS: single
 10 (D) TOPOLOGY: linear
- (ii) MOLECULE TYPE: other nucleic acid

15

(xi) SEQUENCE DESCRIPTION: SEQ ID NO: 2:

20

AAACTTTTAC AGGGATGGAG AACG

24

(2) INFORMATION FOR SEQ ID NO: 3:

- 25 (i) SEQUENCE CHARACTERISTICS:
 (A) LENGTH: 30 base pairs
 (B) TYPE: nucleic acid
 (C) STRANDEDNESS: single
 30 (D) TOPOLOGY: linear
- (ii) MOLECULE TYPE: other nucleic acid

35

(xi) SEQUENCE DESCRIPTION: SEQ ID NO: 3:

40

AGAGGATTAT CTATGCAAAT CCTTGTAACC

30

(2) INFORMATION FOR SEQ ID NO: 4:

45

- (i) SEQUENCE CHARACTERISTICS:
 (A) LENGTH: 30 base pairs
 (B) TYPE: nucleic acid
 50 (C) STRANDEDNESS: single
 (D) TOPOLOGY: linear

- (ii) MOLECULE TYPE: other nucleic acid

55

EP 0 928 832 A2

(xi) SEQUENCE DESCRIPTION: SEQ ID NO: 4:

5

TCAACTTCAC TATCAAAAGT CATCATCTAG

30

(2) INFORMATION FOR SEQ ID NO: 5:

10

(i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 30 base pairs

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

15

(ii) MOLECULE TYPE: other nucleic acid

20

(xi) SEQUENCE DESCRIPTION: SEQ ID NO: 5:

25

TCTTGGGATT CAATACTTT GCAACAGTCA

30

(2) INFORMATION FOR SEQ ID NO: 6:

30

(i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 25 base pairs

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

35

(ii) MOLECULE TYPE: other nucleic acid

40

(xi) SEQUENCE DESCRIPTION: SEQ ID NO: 6:

45

GAATTCCCAA ACTTTTAGAG ACATC

25

(2) INFORMATION FOR SEQ ID NO: 7:

50

(i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 40 base pairs

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

55

(ii) MOLECULE TYPE: other nucleic acid

5

(xi) SEQUENCE DESCRIPTION: SEQ ID NO: 7:

10

TACTAAAAGT GACTCTCTAA TTTCTATTT TTGGAATTA

40

(2) INFORMATION FOR SEQ ID NO: 8:

15

(i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 30 base pairs

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

20

(ii) MOLECULE TYPE: other nucleic acid

25

(xi) SEQUENCE DESCRIPTION: SEQ ID NO: 8:

30

AGTTTGCAGA GAAAGACAAT ATAGTTCTCT

30

(2) INFORMATION FOR SEQ ID NO: 9:

35

(i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 29 base pairs

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

40

(ii) MOLECULE TYPE: other nucleic acid

45

(xi) SEQUENCE DESCRIPTION: SEQ ID NO: 9:

50

TAATCTCTAC CAAATCTGGA TACTATACC

29

(2) INFORMATION FOR SEQ ID NO: 10:

55

(i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 30 base pairs

EP 0 928 832 A2

(B) TYPE: nucleic acid
(C) STRANDEDNESS: single
(D) TOPOLOGY: linear

5

(ii) MOLECULE TYPE: other nucleic acid

10

(xi) SEQUENCE DESCRIPTION: SEQ ID NO: 10:

15

TGATCACTCC ACTGTTTCATA GGGATCCATC

30

(2) INFORMATION FOR SEQ ID NO: 11:

20

(i) SEQUENCE CHARACTERISTICS:
(A) LENGTH: 28 base pairs
(B) TYPE: nucleic acid
(C) STRANDEDNESS: single
(D) TOPOLOGY: linear

25

(ii) MOLECULE TYPE: other nucleic acid

30

(xi) SEQUENCE DESCRIPTION: SEQ ID NO: 11:

35

AATTGAGAG AACTTGATGG TAAGTACA

28

(2) INFORMATION FOR SEQ ID NO: 12:

40

(i) SEQUENCE CHARACTERISTICS:
(A) LENGTH: 30 base pairs
(B) TYPE: nucleic acid
(C) STRANDEDNESS: single
(D) TOPOLOGY: linear

45

(ii) MOLECULE TYPE: other nucleic acid

50

(xi) SEQUENCE DESCRIPTION: SEQ ID NO: 12:

55

GATCTATAT TCATCATAGG AAACACCATT

30

EP 0 928 832 A2

(2) INFORMATION FOR SEQ ID NO: 13:

(i) SEQUENCE CHARACTERISTICS:

- (A) LENGTH: 30 base pairs
- (B) TYPE: nucleic acid
- (C) STRANDEDNESS: single
- (D) TOPOLOGY: linear

(ii) MOLECULE TYPE: other nucleic acid

(xi) SEQUENCE DESCRIPTION: SEQ ID NO: 13:

CCAGACTTCA CTTCTAATGA TGATTATGGG

(2) INFORMATION FOR SEQ ID NO: 14:

(i) SEQUENCE CHARACTERISTICS:

- (A) LENGTH: 27 base pairs
- (B) TYPE: nucleic acid
- (C) STRANDEDNESS: single
- (D) TOPOLOGY: linear

(ii) MOLECULE TYPE: other nucleic acid

(xi) SEQUENCE DESCRIPTION: SEQ ID NO: 14:

ACATTTTCCTT TCAGGGTGTC TGAATAA

(2) INFORMATION FOR SEQ ID NO: 15:

(i) SEQUENCE CHARACTERISTICS:

- (A) LENGTH: 31 base pairs
- (B) TYPE: nucleic acid
- (C) STRANDEDNESS: single
- (D) TOPOLOGY: linear

(ii) MOLECULE TYPE: other nucleic acid

EP 0 928 832 A2

(xi) SEQUENCE DESCRIPTION: SEQ ID NO: 15:

TTGTGGATCA AATTTCAGTT GACTTGTCAT C

31

5

(2) INFORMATION FOR SEQ ID NO: 16:

(i) SEQUENCE CHARACTERISTICS:

10

(A) LENGTH: 30 base pairs

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

15

(ii) MOLECULE TYPE: other nucleic acid

20

(xi) SEQUENCE DESCRIPTION: SEQ ID NO: 16:

GTATCTATAT TCATCATAGG AAACACCACA

30

25

(2) INFORMATION FOR SEQ ID NO: 17:

(i) SEQUENCE CHARACTERISTICS:

30

(A) LENGTH: 30 base pairs

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

35

(ii) MOLECULE TYPE: other nucleic acid

40

(xi) SEQUENCE DESCRIPTION: SEQ ID NO: 17:

GACTTCACTT CTAATGATGA TTATGGGAGA

30

45

(2) INFORMATION FOR SEQ ID NO: 18:

(i) SEQUENCE CHARACTERISTICS:

50

(A) LENGTH: 30 base pairs

(B) TYPE: nucleic acid

(C) STRANDEDNESS: single

(D) TOPOLOGY: linear

55

(ii) MOLECULE TYPE: other nucleic acid

EP 0 928 832 A2

(xi) SEQUENCE DESCRIPTION: SEQ ID NO: 18:

5 TGCCATGGGG CCTGTGCAAG GAAGTATTGA 30

(2) INFORMATION FOR SEQ ID NO: 19:

10 (i) SEQUENCE CHARACTERISTICS:
(A) LENGTH: 30 base pairs
(B) TYPE: nucleic acid
(C) STRANDEDNESS: single
15 (D) TOPOLOGY: linear

(ii) MOLECULE TYPE: other nucleic acid

20

(xi) SEQUENCE DESCRIPTION: SEQ ID NO: 19:

25 AGCCTATGCC TAGATAAATC GCGATAGACT 30

(2) INFORMATION FOR SEQ ID NO: 20:

30 (i) SEQUENCE CHARACTERISTICS:
(A) LENGTH: 37 base pairs
(B) TYPE: nucleic acid
(C) STRANDEDNESS: single
35 (D) TOPOLOGY: linear

(ii) MOLECULE TYPE: other nucleic acid

40

(xi) SEQUENCE DESCRIPTION: SEQ ID NO: 20:

45 GTTTCACATA GTGTATGACC CTCTATATAC ACTCATT 37

(2) INFORMATION FOR SEQ ID NO: 21:

50

(i) SEQUENCE CHARACTERISTICS:
(A) LENGTH: 30 base pairs
(B) TYPE: nucleic acid

55

EP 0 928 832 A2

(C) STRANDEDNESS: single
(D) TOPOLOGY: linear

5 (ii) MOLECULE TYPE: other nucleic acid

10 (xi) SEQUENCE DESCRIPTION: SEQ ID NO: 21:

CCTATGCACT AATCAAAGGA ATCATCCTGT

30

15 (2) INFORMATION FOR SEQ ID NO: 22:

(i) SEQUENCE CHARACTERISTICS:
(A) LENGTH: 30 base pairs
(B) TYPE: nucleic acid
(C) STRANDEDNESS: single
(D) TOPOLOGY: linear

20

25 (ii) MOLECULE TYPE: other nucleic acid

30

(xi) SEQUENCE DESCRIPTION: SEQ ID NO: 22:

TTTGTTTATT GCTCCAAGAG AGTCATACCA

30

35 (2) INFORMATION FOR SEQ ID NO: 23:

(i) SEQUENCE CHARACTERISTICS:
(A) LENGTH: 25 base pairs
(B) TYPE: nucleic acid
(C) STRANDEDNESS: single
(D) TOPOLOGY: linear

40

45 (ii) MOLECULE TYPE: other nucleic acid

50

(xi) SEQUENCE DESCRIPTION: SEQ ID NO: 23:

GCTAAAGAAA TTCTTGCTCG TTGTT

25

55

(2) INFORMATION FOR SEQ ID NO: 24:

(i) SEQUENCE CHARACTERISTICS:

- (A) LENGTH: 60 base pairs
- (B) TYPE: nucleic acid
- (C) STRANDEDNESS: single
- (D) TOPOLOGY: linear

(ii) MOLECULE TYPE: other nucleic acid

(xi) SEQUENCE DESCRIPTION: SEQ ID NO: 24:

GACTGACTGA CTGACTGACT CTGACTGACT TATTCACCTT GCTAAAGAAA TTCTTGCTGA 60

(2) INFORMATION FOR SEQ ID NO: 25:

(i) SEQUENCE CHARACTERISTICS:

- (A) LENGTH: 28 base pairs
- (B) TYPE: nucleic acid
- (C) STRANDEDNESS: single
- (D) TOPOLOGY: linear

(ii) MOLECULE TYPE: other nucleic acid

(xi) SEQUENCE DESCRIPTION: SEQ ID NO: 25:

TAAAATTGGA GCAATGTTGT TTTTGACC 28

(2) INFORMATION FOR SEQ ID NO: 26:

(i) SEQUENCE CHARACTERISTICS:

- (A) LENGTH: 30 base pairs
- (B) TYPE: nucleic acid
- (C) STRANDEDNESS: single
- (D) TOPOLOGY: linear

(ii) MOLECULE TYPE: other nucleic acid

(xi) SEQUENCE DESCRIPTION: SEQ ID NO: 26:

5 TATTTTATT TCAGATGCGA TCTGTGAGTT 30

(2) INFORMATION FOR SEQ ID NO: 27:

10 (i) SEQUENCE CHARACTERISTICS:
 (A) LENGTH: 29 base pairs
 (B) TYPE: nucleic acid
 (C) STRANDEDNESS: single
 15 (D) TOPOLOGY: linear

(ii) MOLECULE TYPE: other nucleic acid

20

(xi) SEQUENCE DESCRIPTION: SEQ ID NO: 27:

25 TTTTGCTGTG AGATCTTTGA CAGTCATTT 29

30 Claims

1. A method for detecting the presence or absence of twelve mutations in the cystic fibrosis transmembrane conductor regulator (CFTR) gene, which method comprises contacting sample genomic DNA from an individual in two separate reaction vessels with allele specific primer sets for (A) 1717-1 G>A, G542X, W1282X, N1303K, ΔF508(M), 3849+10kb C>T mutations and (B) the 621+1 G>T, R553X, G551D, R117H, R1162X and R334W mutations respectively, in the presence of appropriate nucleotide triphosphates and an agent for polymerisation, such that each diagnostic primer is extended only when the relevant mutation is present in the sample; and detecting the presence or absence of CFTR gene alleles by reference to the presence or absence of diagnostic primer extension product (s).
2. A method as claimed in claim 1 and wherein one or more diagnostic primers is used with one or more amplification primers in one or more cycles of PCR amplification.
3. A set of allele specific primers for each of the following alleles of the CFTR gene: 1717-1 G>A, G542X, W1282X, N1303K, ΔF508(M), and 3849+10kb C>T mutations.
4. A set of primers as claimed in claim 3 and comprising the following diagnostic primer sequences:

50 TCTTGGGATTCAATAACTTTGCAACAGTCA
 TACTAAAAGTGACTCTCTAATTTTCTATTTTGGTAATTA
 AGTTTGACAGAGAAAGACAATATAGTTCTCT
 TGATCACTCCACTGTTTCATAGGGATCCATC
 55 GTATCTATATTCATCATAGGAAACACCATT
 ACATTTTCCTTTCAGGGTGTCTGACTAA

EP 0 928 832 A2

5. A set of allele specific primers for each of the following alleles of the CFTR gene: 621+1 G>T, R553X, G551D, R117H, R1162X and R334W mutations.

6. A set of primers as claimed in claim 5 and comprising the following diagnostic primer sequences:

5

10

15

GTATCTATATTCATCATAGGAAACACCACA
TGCCATGGGGCCTGTGCAAGGAAGTATTGA
AGCCTATGCCTAGATAAATCGCGATAGACT
CCTATGCACTAATCAAAGGAATCATCCTGT
GCTAAAGAAATTCTTGCTCGTTGTT
GACTGACTGACTGACTGACTCTGACTGACTTATTCA
CCTTGCTAAAGAAATTCTTGCTGA
TATTTTATTTTCAGATGCGATCTGTGAGTT

20

7. A set of primers comprising the following diagnostic primer and amplification primer sequences:

25

30

35

TCTTGGGATTCAATAACTTTGCAACAGTCA
GAATTCCTAACTTTAGAGACATC
TACTAAAAGTGACTCTCTAATTTTCTATTTTGGTAATTA
AGTTTGAGAGAAAGACAATATAGTTCTCT
TAATCTCTACCAAATCTGGATACTATACC
TGATCACTCCACTGTTTCATAGGGATCCATC
AATTTGAGAGAACTTGATGGTAAGTACA
GTATCTATATTCATCATAGGAAACACCATT
CCAGACTTCACTTCTAATGATGATTATGGG
ACATTTCTTTTCAGGGTGTCTGACTAA
TTGTGGATCAAATTTAGTTGACTTGTCTATC

40

8. A set of primers comprising the following diagnostic primer and amplification primer sequences:

45

50

55

GTATCTATATTCATCATAGGAAACACCACA
GACTTCACTTCTAATGATGATTATGGGAGA
TGCCATGGGGCCTGTGCAAGGAAGTATTGA
AGCCTATGCCTAGATAAATCGCGATAGACT
GTTTCACATAGTGTATGACCTCTATATACACTCATT
CCTATGCACTAATCAAAGGAATCATCCTGT
TTTGTTTATTGCTCCAAGAGAGTCATACCA
GCTAAAGAAATTCTTGCTCGTTGTT
GACTGACTGACTGACTGACTCTGACTGACTTATTCA

EP 0 928 832 A2

CCTTGCTAAAGAAATTCTTGCTGA
TAAAATTGGAGCAATGTTGTTTTGACC
TATTTTATTTTATTCAGATGCGATCTGTGAGTT
TTTTGCTGTGAGATCTTTGACAGTCATT

5

10

9. A set of primers as claimed in any one of the previous claims and comprising one or more of the following control primers:

GAGCACAGTACGAAAAACCACT
AAACTTTTACAGGGATGGAGAACG
AGAGGATTATCTATGCAAATCCTTGTAACC
TCAACTTCACTATCAAAAGTCATCATCTAG

15

20

10. A diagnostic kit for detecting the presence or absence of twelve mutations in the cystic fibrosis transmembrane conductor regulator (CFTR) gene which comprises sets of primers as claimed in any of the previous claims.

25

30

35

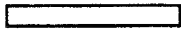


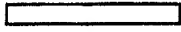
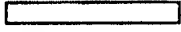


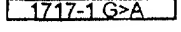
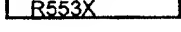

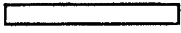
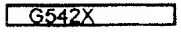

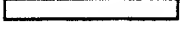
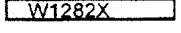
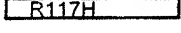
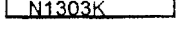
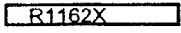
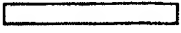
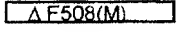

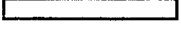
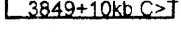
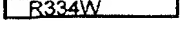


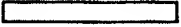
40

45

50

55

FIG 1

A-tube	B-tube	Marker
		 (550)
 (487)	 (487)	 (500)
		 (450)
		 (400)
	 (383)	
 (329)	 (325)	 (350)
		 (300)
 (279)	 (285)	
		 (250)
 (240)	 (243)	
 (212)	 (200)	 (200)
 (160)	 (160)	
		 (150)
 (130)	 (140)	
 (97)	 (97)	 (100)



DOCUMENTS CONSIDERED TO BE RELEVANT			
Category	Citation of document with indication, where appropriate, of relevant passages	Relevant to claim	CLASSIFICATION OF THE APPLICATION (Int.Cl.7)
Y	WO 92 01066 A (GENETYPE AG) 23 January 1992 (1992-01-23) * page 3, line 6 - line 12 * * page 12, line 10 - line 12 * * page 18, line 9 - line 10 * * page 49, line 32 - line 37 * * page 50, line 1 - line 3 * * page 56, line 3 - line 10 * * page 57, line 12 - line 18 * * page 69, line 2 - line 14 * * page 69, line 30 - line 34 * ---	1-3	C12Q1/68
Y	EP 0 928 832 A (ZENECA LTD) 14 July 1999 (1999-07-14) * page 2, line 31 - line 32 * * page 2, line 53 - line 54 * * page 2, line 51 * ---	1-3	
Y	JP 10 084999 T (HITACHI ELECTRON ENG CO LTD) 7 April 1998 (1998-04-07) * abstract * * claim 1 * * paragraphs '0028!', '0029!' * -----	3	
			TECHNICAL FIELDS SEARCHED (Int.Cl.7)
			C12Q
The present search report has been drawn up for all claims			
Place of search MUNICH		Date of completion of the search 5 July 2002	Examiner Helliot, B
CATEGORY OF CITED DOCUMENTS X : particularly relevant if taken alone Y : particularly relevant if combined with another document of the same category A : technological background O : non-written disclosure P : intermediate document T : theory or principle underlying the invention E : earlier patent document, but published on, or after the filing date D : document cited in the application L : document cited for other reasons & : member of the same patent family, corresponding document			

**ANNEX TO THE EUROPEAN SEARCH REPORT
ON EUROPEAN PATENT APPLICATION NO.**

EP 01 30 6444

This annex lists the patent family members relating to the patent documents cited in the above-mentioned European search report. The members are as contained in the European Patent Office EDP file on
The European Patent Office is in no way liable for these particulars which are merely given for the purpose of information.

05-07-2002

Patent document cited in search report		Publication date	Patent family member(s)		Publication date
WO 9201066	A	23-01-1992	AT	185377 T	15-10-1999
			AU	647806 B2	31-03-1994
			AU	8099491 A	04-02-1992
			WO	9201066 A1	23-01-1992
			CA	2087042 A1	12-01-1992
			DE	69131691 D1	11-11-1999
			DE	69131691 T2	31-05-2000
			DK	570371 T3	17-04-2000
			EP	0570371 A1	24-11-1993
			IE	912426 A1	15-01-1992
			IL	98793 A	31-10-1995
			NZ	238926 A	23-12-1993
			US	5851762 A	22-12-1998
			ZA	9105422 A	13-04-1993
EP 0928832	A	14-07-1999	EP	0928832 A2	14-07-1999
JP 10084999	T		NONE		